MANUAL

Salente®

SlimFit

CZ

SK

EN

HU

RO

DE

IT

Obsah / Content / Tartalom / Conținut / Inhaltsübersicht / Indice

Česky 4

Slovensky 16

English 28

Magyar 40

Română 52

Deutsch 64

Italiano 77

- CZ
- SK
- EN
- HU
- RO
- DE
- IT

Obsah

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Table of Contents

Home..... 28

THE MOST UP-TO-DATE VERSION OF THE USER MANUAL..... 28

SAFETY WARNINGS AND INSTRUCTIONS 28

PACKAGING CONTENT 30

DESCRIPTION OF PARTS 31

SPECIFICATIONS 32

PRINCIPLE OF MEASUREMENT 32

SERVICE 33

LINKING THE SCALE TO THE APP 35

ERROR MESSAGES..... 37

BATTERY CHANGE 37

CLEANING AND MAINTENANCE 37

- Before putting this product into operation, please read this entire user manual carefully through to the end, even if you are familiar with the use of a similar type of product in the past.
- Keep this user manual for future reference - it is an integral part of the product and may contain important instructions on commissioning, operation and care of the product.
- Make sure that other people who will use this product read this user manual.
- If you give, gift or sell this product to others, please give them this user manual together with this product and make them aware of the need to read it.
- During the warranty period, we recommend that you keep the original packaging of the product - this way the product will be protected in the best possible way in case of transport.

THE MOST UP-TO-DATE VERSION OF THE USER MANUAL

Due to Salente's ongoing product innovations and changes in legislation, you can always find the most up-to-date version of this user manual on the website www.salente.eu or by simply scanning the QR code shown below with your mobile phone:

**SAFETY WARNINGS AND INSTRUCTIONS**

- Use the product only as specified in this user manual.
- This product may be used by children aged 8 years and older and persons with reduced physical, sensory or mental abilities or lack of experience and knowledge, provided they are supervised or have been instructed in the safe use of the appliance and understand the potential hazards.
- Cleaning and maintenance by the user must not be carried out by unsupervised children.

- Unpack the product carefully and be careful not to throw away any part of the packaging material before you have found all the components of the product.
- Remove all packaging and marketing materials from the product before first use.
- Exercise extra caution if the product will be used around children. Always place the product out of their reach. Children should be supervised to ensure that they do not play with the product.
- Do not allow the product to be used as a toy.
- Always install the product on a flat, hard, dry and stable place.
- Never place the product on inclined surfaces.
- Do not use the product other than as described in this user manual.
- The product is designed for use in homes, offices and similar areas. Do not use it for commercial, medical or other purposes.
- The product is not intended for continuous operation.
- This product is designed for indoor use, do not use it outdoors.
- Do not use the product within the range of a strong electromagnetic field (e.g. microwave oven).
- Never use accessories that are not supplied with or intended for this product.
- Never repair the product yourself or make any modifications to it! Entrust all repairs and adjustments to this product to a professional company/service. Tampering with the product during the warranty period puts you at risk of losing warranty benefits.
- Do not shake the product and take care not to drop it on the ground.
- Never use the product if it is not working properly, if it has fallen on the ground and become damaged, or if it has fallen into water. In this case, have it repaired by a professional company/service.
- Use and store the product away from explosive gases, flammable, combustible and volatile substances (solvents, paints, adhesives, etc.).
- Do not expose the product to extreme temperatures, direct sunlight, excessive humidity (e.g. sauna) or place it in excessively dusty environments.
- Do not place the product near heaters, open flames or other appliances or equipment that are sources of heat.
- Do not use the product for purposes other than those for which it is intended.
- Do not spray or rinse the product with water or other liquid, or immerse it in water or other liquid, even partially.
- Before stepping onto the weighing surface, make sure that the surface is not wet or damp and that your feet are not wet. Otherwise, there is a risk of slipping and injury.
- Do not jump on the scale and do not exceed its maximum weighing capacity. Otherwise it could be damaged.


- Children under 10 years of age and persons over 100 years of age may only use this product in the weighing mode.
- The manufacturer is not liable for any damage caused by improper use of the appliance and its accessories.
- If you will not be using the product for a longer period of time (more than 3 months), remove the batteries from the product, otherwise there is a risk of them leaking into the battery compartment and irreparably damaging the product.
- Keep the batteries out of the reach of children. In their presence, they can pose a danger. If the child swallows the battery, seek medical attention immediately.
- Do not use the product to store objects and ensure that no other objects are placed on the product when storing it.
- Keep the product clean. When cleaning, follow the instructions in the Cleaning and Maintenance chapter.
- Cleaning and maintenance by the user must not be carried out by children unless they are at least 8 years old and supervised. Children under 8 years of age must stay away from the product.

PACKAGING CONTENT


- Smart personal scale
- AAA micro-core alkaline batteries (3 pcs)
- Operating Instructions

DESCRIPTION OF PARTS

Weight:

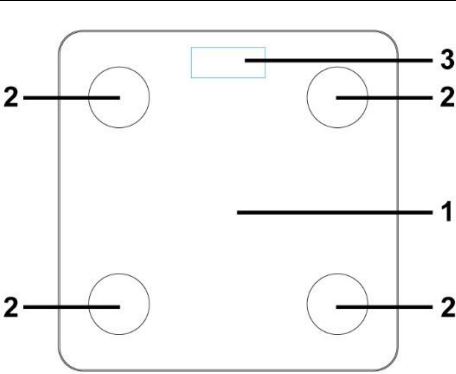 Diagram showing the top view of the scale platform. Callout 1 points to the central weighing area. Callout 2 points to four circular foot sensors (electrodes) located at the corners. Callout 3 points to a small rectangular digital display at the top center. Callout 4 points to a button on the right side. Callout 5 points to a rectangular battery compartment cover on the right side. Callout 6 points to four circular non-slip feet with integrated weighing sensors at the corners.	4) Button for switching units (kg/lb/hr) 5) Battery compartment cover 6) Non-slip feet with integrated weighing sensors
1) Weighing area 2) Foot sensors (electrodes) 3) Digital display	

Display:

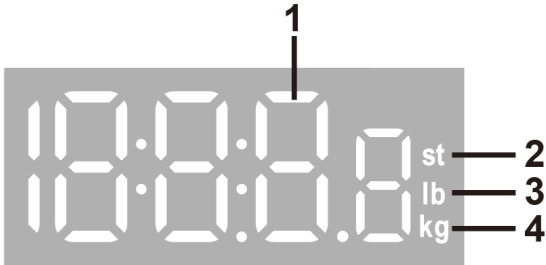 Diagram showing the digital display. Callout 1 points to the measured weight indicator (00.00). Callout 2 points to the unit of mass in the imperial system (st). Callout 3 points to the unit of mass in the imperial system (lb). Callout 4 points to the unit of weight in the metric system (kg).	1) Measured weight indicator 2) Indication of the unit of mass in the imperial system (st)
	3) Indication of the unit of mass in the imperial system (lb) 4) Indication of the unit of weight in the metric system (kg)

SPECIFICATIONS

- Max. weighing capacity 180 kg
- Resolution 100 g
- Measurement and analysis of 15 parameters (weight, BMI index, % fat, muscle, water, etc.)
- Data history storage in the MovingLife mobile app with support for 23 languages (CZ, SK, EN, DE, HU, PL, RO, etc.)
- Support iOS (Apple) from version 8.0 and Android from version 4.3
- Up to 8 mobile devices can be paired
- Connectivity with Apple Health and Fitbit
- Illuminated MagicVision display with extra large digits (2.1 cm high)
- 4 high-precision "StrainGauge" weighing sensors
- Weighing surface made of 5 mm thick safety glass
- Double AutoOn function (automatic switching on of the scale when entering the weighing area and immediate start of weighing)
- Choice of metric or imperial units (kg, lb, st)
- BatterySaver (automatic shutdown of the scale after 30 seconds of inactivity)
- Warning indication of exceeding the maximum weighing capacity of the scale
- Warning indication when battery needs to be replaced
- Body fat measurement range: 4-75 %
- Accuracy of body fat measurement: 0,1 %
- The range of water measurements in the body: 15-66,5 %
- Accuracy of water measurement in the body: 0,1 %
- Muscle tissue measurement range: 23-91 %
- Accuracy of muscle tissue measurement: 0,1 %
- Bone tissue measurement range: 1-6 %
- Accuracy of bone tissue measurements: 0,1 %
- Accuracy of the calculation of the recommended daily calorie intake (BMR): 1 kcal
- Power supply: 3x 1.5 V AAA micro alkaline batteries (included)
- Certification CE, RoHS
- Dimensions (w × h × h): 28 × 28 × 2.2 cm
- Weight (without batteries): 1108 g

PRINCIPLE OF MEASUREMENT

This scale uses bioimpedance analysis (BIA) to measure and analyse body parameters (percentage of bone, fat and muscle tissue, water in the body, etc.) based on the principle of a weak electrical impulse sent through the scale's foot sensors (electrodes) into the human body. The software embedded in the scale evaluates the

measurements based on the difference in electrical impedance of the individual tissues and their relative contribution to the body.

Muscle contains about 90% water and adipose tissue 45% water. Therefore, the percentage sum can exceed 100% and for this reason the measurements cannot be added because both muscle and fat contain water.

SERVICE

Important notices:

- If the scale is placed on an uneven or soft surface (e.g. carpet), inaccuracies in weighing may occur (the scale may display higher or lower readings).
- The use of body parameter analysis mode (i.e. barefoot weighing on sensors) is not recommended:
 - pregnant women
 - people with an implanted pacemaker or other medical implant
 - people with symptoms of oedema (swelling)
 - persons undergoing dialysis treatment

Preparation for weighing

Warning:

The batteries are already inserted in the battery compartment at the factory. To prevent the batteries from discharging during transport from the factory to the customer, the batteries are protected by an insulating strip with the inscription "REMOVE BEFORE USE", which must be removed from the battery compartment before the first use by pulling it.

1. After removing the insulating tape, place the scale on a flat, dry and hard surface.
2. The symbol "8888" will briefly appear on the display to indicate the stabilization of the weighing sensors into weighing mode.
3. After the display shows "0.0 kg", the weighing can start.
4. If necessary, switch the units using the switch at the bottom of the scale (the scale is factory set to metric units or weighing in kilograms).
5. The measured weight remains on the display for 10 seconds even after stepping off the scale.

Notes:

- The scale is factory preset to weigh in kilograms (the "kg" symbol appears next to the measured weight on the display).

- The scale is equipped with the "Double Auto On" function, where it automatically switches itself on as soon as you step on it and starts weighing the weight immediately.

Weight weighing mode

In this mode, the scale weighs only the weight.

1. Step on the scale with your feet so that the entire surface of your feet touch the weighing surface. Do not stand on your toes or heels only, do not lean and stand still until the final weight is displayed.
2. The weighing progress is indicated by a triple flash of the value on the display.
3. The resulting weight remains displayed for 10 seconds. Then the display automatically switches off and after a further 20 seconds of inactivity the whole scale automatically switches off.

Warning:

- We recommend weighing on the same pad/floor and at the same time of day, ideally first thing in the morning on an empty stomach - this will give comparable results.

Weighing mode with body parameters measurement

In this mode, the scale weighs the weight and evaluates the body parameters at the same time, and it is necessary to stand on the scale with bare feet during weighing so that the toe touches the upper and the heel touches the lower foot sensor.

1. Step on the scale with bare feet so that the entire surface of the feet touches the weighing surface. Do not stand on your toes or heels only, do not lean and stand still until the final weight is displayed.
2. The weighing progress is indicated by the successive flashing of four "o" characters below the weight value on the display.
3. The resulting weight remains displayed for 10 seconds. Then the display automatically switches off and after a further 20 seconds of inactivity the whole scale automatically switches off.
4. The values are then transferred to the application where they are analyzed and displayed.

Warning:

- Measurement of body parameters (except weight) is not recommended immediately after waking up, as the body is dehydrated after sleep. Similarly,

we do not recommend measurements immediately after sports performance or during drastic weight loss diets.

- Wet, dirty feet or feet covered in cream can distort the results of body measurements or even make them impossible.

LINKING THE SCALE TO THE APP

This scale allows you to transfer, analyse and display all the measured values in the MovingLife mobile app. The values are transferred to the app via Bluetooth.

Sharing Libra with multiple users


1. Open App > Me > Scale with WiFi > Share QR code
2. In the next step, the second user scans the QR code to make the connection.
3. Support for up to eight users simultaneously

How to pair your mobile phone with a smart scale

Turn on Bluetooth and connect your device to the internet (using 4G or Wi-Fi), launch the MovingLife app, and enter Account→My Device→Add Device. Then select the appropriate device and pair them.

Account:

1. Local
 - a. No registration is required, all data will only remain on your device. If you uninstall the app or start using another device, you will lose data.
2. Cloud
 - a. By switching the switch to the on position, you can register to the Cloud and the data from the app will remain there even if you replace your device.

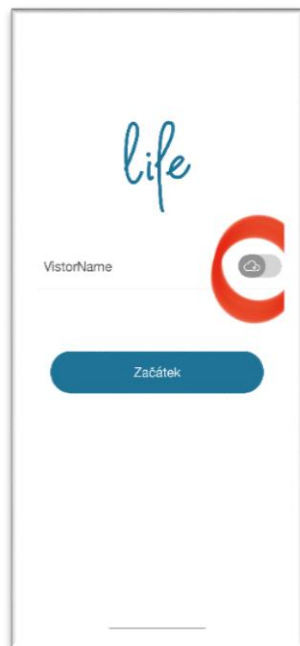For the app to work properly and count health functions, you must fill in the basic user information that the app asks you for during initial setup. The app can be linked to Apple Health or Google Fit, and synchronisation can be enabled in the app: "Me"→"Settings" and select Google Fit or Apple Health depending on your device

Pairing failed

Please check the following items:

- 1) Verify that your mobile phone's operating system and Bluetooth version meet the conditions.
- 2) Stand on the scale and check if the display has lit up. If not, change the batteries in time.
- 3) Make sure your mobile phone has Bluetooth capability.
- 4) If all the above items are OK, restart the app or restart your phone.

Warning:

The app requires Bluetooth version 4.2 support and, for mobile phones, Apple iOS 12.5 and above or Android 9.0 and above.

Warning:

- If the scale is placed on an uneven or soft surface (e.g. carpet), inaccuracies in weighing may occur (the scale may display higher or lower readings).
- Wet, dirty feet or feet covered in cream can distort the results of body measurements or even make them impossible.
- We do not recommend measuring body parameters immediately after waking up, as the body is dehydrated after sleep. Similarly, we do not recommend measurements immediately after sports performance or during drastic weight loss diets.
- Variations in measured physical parameters may be greater in children under 10 years of age, persons over 70 years of age and professional athletes.
- The measured values of the body parameters are only indicative and no medical conclusions can be drawn from them. If the results are significantly different from normal, we recommend consulting a nutritionist or doctor.
- The weight unit you set remains stored in the scale's internal memory even after the scale is switched off and on again. The exception is when removing or replacing the battery, when the scale is switched back on, it displays the unit for weighing the weight in grams (the g symbol appears on the display).

Options for switching off the scale

If the scale remains idle for more than 10 seconds, the display will automatically switch off. After 30 seconds of inactivity, the entire scale automatically shuts down.

ERROR MESSAGES

Lo... low battery indication. Replace the batteries with new ones.

Err... the maximum weighing capacity (180 kg) has been exceeded. Get off the scale immediately, otherwise the weighing sensors or the entire product may be irreversibly damaged.

-Err... the scale sensors were loaded before switching on. Wait about 2 seconds for the scale to reset itself.

C... there was a weighing error. Repeat the weighing.

BATTERY CHANGE

1. Open the battery compartment cover at the bottom of the scale.
2. Remove the old batteries and insert the new batteries into the battery compartment. Make sure to observe the correct polarity, which is marked on the bottom of the battery compartment.
3. Close the battery compartment by snapping its cover back on.

Warning:

- Do not use rechargeable batteries to power the scale, their voltage may not be sufficient to power the scale.
- Do not combine old and new batteries due to their different residual capacities. A low battery error indication could be displayed prematurely and all batteries may need to be replaced.

CLEANING AND MAINTENANCE

- Never clean the scale under running water, rinse it or immerse it in water.
- Use a soft cloth dampened with lukewarm water to clean the outside of the scale. Do not use harsh cleaners, thinners or solvents or the scale surface may be damaged.

— **Disposal:** If the information leaflet, packaging or product label contains a crossed-out bin symbol, this means that - in the European Union - all electrical and electronic products, batteries and accumulators marked in this way must be placed in separate waste bins at the end of their life. Do not dispose of these products in unsorted municipal waste.

CE DECLARATION OF CONFORMITY

Abacus Electric, s.r.o. hereby declares that this product meets the requirements of the standards and regulations relevant to this type of device, and also declares that the product complies with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <http://ftp.salente.cz/ce>

Importer / Producer

Abacus Electric, s.r.o.

Planá 2, 370 01, Czech Republic

E-mail: helpdesk@salente.eu

Made in P.R.C.

Copyright © Abacus Electric, s.r.o.

www.salente.eu

All rights reserved.

Design and appearance are subject to change without notice.

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