MANUAL

Salente®

BodyFit

CZ

SK

EN

HU

RO

DE

SI

Česky 4

Slovensky 14

English 24

Magyar 34

Română 44

Deutsch 54

Slovinčina 65

- CZ
- SK
- EN
- HU
- RO
- DE
- SI

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Table of Contents

INTRODUCTION 24

THE MOST UP-TO-DATE VERSION OF THE USER MANUAL 24

SAFETY WARNINGS AND INSTRUCTIONS 24

PACKAGING CONTENT 26

DESCRIPTION OF PARTS 26

INFORMATION AND ERROR MESSAGES..... 28

SPECIFICATIONS 28

PURPOSE OF USE..... 29

PRINCIPLE OF MEASUREMENT 29

SERVICE 29

MOVINGLIFE APP 31

CLEANING AND MAINTENANCE 32

EN

INTRODUCTION

WELCOME TO THE SALENTE FAMILY

- Before putting this product into operation, please read this entire user manual carefully through to the end, even if you are familiar with the use of a similar type of product in the past.
- Keep this user manual for future reference - it is an integral part of the product and may contain important instructions on commissioning, operation and care of the product.
- Make sure that other people who will use this product read this user manual.
- If you give, gift or sell this product to others, please give them this user manual together with this product and remind them to read it.
- During the warranty period, we recommend that you keep the original packaging of the product - this way the product will be protected in the best possible way in case of transport.

THE MOST UP-TO-DATE VERSION OF THE USER MANUAL

Due to Salente's ongoing product innovations and changes in legislation, you can always find the most up-to-date version of this user manual on the website www.salente.eu or by simply scanning the QR code shown below with your mobile phone:

SAFETY WARNINGS AND INSTRUCTIONS

- **This appliance may be used by children aged 8 years or older and persons with reduced physical, sensory or mental abilities or lack of experience and knowledge, provided they are supervised or have been instructed in the safe use of the appliance and understand the potential hazards.**
- **Children must not play with the appliance.**
- **Do not expose the appliance to dripping or splashing water. Do not spray or rinse it with water or other liquid. Do not immerse it in water or any other liquid, even partially.**

- The appliance is designed for indoor use, do not use it outdoors.**
- When cleaning, follow the instructions in the Cleaning and Maintenance chapter.**
- Cleaning and maintenance by the user must not be carried out by children unless they are over 8 years of age and supervised. Children under 8 years of age must stay away from the product.**
- Use the product only as specified in this user manual.
- Unpack the product carefully and be careful not to throw away any part of the packaging material before you have found all the components of the product.
- Remove all packaging and marketing materials from the product before first use.
- Exercise extra caution if the product will be used around children. Always place the product out of their reach. Children should be supervised to ensure that they do not play with the product.
- Do not allow the product to be used as a toy.
- Always install the product on a flat, hard, dry and stable place.
- Never place the product on inclined surfaces.
- Do not use the product in any other way than described in this user manual.
- The product is designed for use in homes, offices and similar areas. Do not use it for commercial, medical or other purposes.
- The product is not intended for continuous operation.
- This product is designed for indoor use, do not use it outdoors.
- Do not use the product within the range of a strong electromagnetic field (e.g. microwave oven).
- Never use accessories that are not supplied with or intended for this product.
- Never repair the product yourself and do not make any modifications to it! Entrust all repairs and adjustments to this product to a professional company/service. Tampering with the product during the warranty period puts you at risk of losing warranty benefits.
- Do not shake the product and take care not to drop it on the ground.
- Never use the product if it is not working properly, if it has fallen on the ground and become damaged, or if it has fallen into water. In this case, have it repaired by a professional company/service.
- Use and store the product away from explosive gases, flammable, combustible and volatile substances (solvents, paints, adhesives, etc.).
- Do not expose the product to extreme temperatures, direct sunlight, excessive humidity (e.g. sauna) or place it in excessively dusty environments.
- Do not place the product near heaters, open flames or other appliances or equipment that are sources of heat.
- Do not use the product for purposes other than those for which it is intended.
- Do not spray or rinse the product with water or other liquid, or immerse it in water or other liquid, even partially.
- Before stepping onto the weighing surface, make sure that the surface is not wet or damp and that your feet are not wet. Otherwise, there is a risk of slipping and injury.
- Do not jump on the scale and do not exceed its maximum weighing capacity. Otherwise it could be damaged.
- Children under 10 years of age and persons over 100 years of age may only use this product in the weighing mode.
- The manufacturer is not liable for any damage caused by improper use of the appliance and its accessories.
- If you will not be using the product for a longer period of time (more than 3 months), remove the batteries from the product, otherwise there is a risk of them leaking into the battery compartment and irreparably damaging the product.

- Keep the batteries out of the reach of children. In their presence, they can pose a danger. If the child swallows the battery, seek medical attention immediately.
- Do not use the product to store objects and ensure that no other objects are placed on the product when storing it.
- Keep the product clean. When cleaning, follow the instructions in the Cleaning and Maintenance chapter.
- Cleaning and maintenance by the user must not be carried out by children unless they are at least 8 years old and supervised. Children under 8 years of age must stay away from the product.

PACKAGING CONTENT


- Personal scale
- AAA batteries (4 pcs, LR03)
- Operating Instructions

DESCRIPTION OF PARTS


Weight:

The diagram illustrates the components of the personal scale. On the left, a top-down view shows the weighing area (1) and the digital display (2). On the right, a bottom-up view shows the battery compartment cover (4) and the four non-slip feet with integrated weighing sensors (5). A switch (3) is also indicated on the bottom view.

1) Weighing area	4) Battery compartment cover
2) Digital display	5) Non-slip feet with integrated weighing sensors
3) Button for switching units (kg/lb/st)	

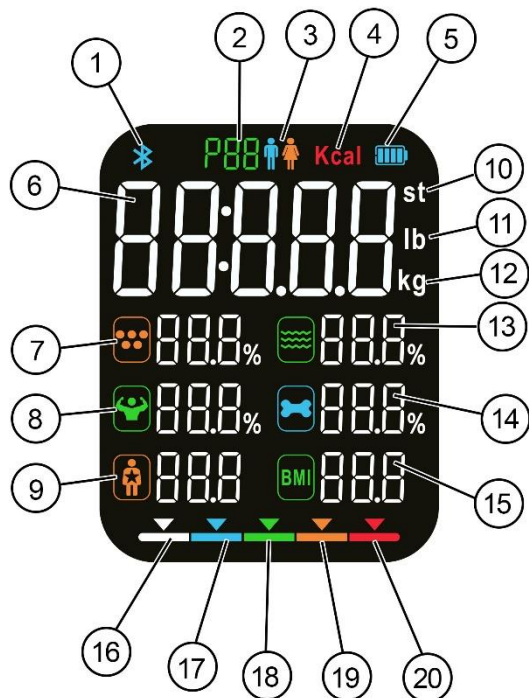- | | |
|--|--|
| 1) Bluetooth active indication | 11) Indication of the unit of mass in the imperial system (lb) |
| 2) User code | 12) Indication of the unit of weight in the metric system (kg) |
| 3) Gender indication (male/female) | 13) Body water percentage indicator |
| 4) Energy unit indication - kilocalorie (Kcal) | 14) An indicator of the percentage of bone tissue in the body |
| 5) Battery status indicator | 15) BMI indicator |
| 6) Numeric display | 16) Graphic indicator of underweight |
| 7) Body fat percentage indicator | 17) Graphical indicator of mild underweight |
| 8) An indicator of the percentage of muscle tissue in the body | 18) Normal weight graphic indicator |
| 9) Metabolic age indicator | 19) Graphical indicator of overweight |
| 10) Indication of the unit of mass in the imperial system (st) | 20) Graphical obesity indicator |

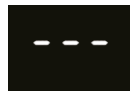	Measurements in progress.
	Low battery indication (replace batteries with new ones).
	Exceeding the maximum weighing capacity (180 kg). Get off the weighing surface immediately, otherwise irreparable damage to the weighing sensors or the entire product may occur.
	Measurement error (get off the scale and repeat the measurement).

SPECIFICATIONS

- Max. weighing capacity 180 kg
- Resolution: 10 g (range 100 g to 10 kg), 50 g (range 10 kg to 180 kg)
- Measurement and analysis of 15 parameters (weight, BMI index, % fat, muscle, water, etc.)
- Data history storage in the MovingLife mobile app with support for 23 languages (CZ, SK, EN, DE, HU, PL, RO, etc.)
- Support iOS (Apple) from version 8.0 and Android from version 4.3
- Up to 8 mobile devices can be paired
- Connectivity with Apple Health and Fitbit
- Backlit MagicVision display with extra large digits (2.5 cm high)
- 4 high-precision "StrainGauge" weighing sensors
- Weighing surface made of 6 mm thick safety glass
- Double AutoOn function (automatic switching on of the scale when entering the weighing area and immediate start of weighing)
- Choice of metric or imperial units (kg, lb, st)
- BatterySaver (automatic shutdown of the scale after 30 seconds of inactivity)
- Warning indication of exceeding the maximum weighing capacity of the scale
- Warning indication when battery needs to be replaced
- Body fat measurement range: 4-75 %
- Accuracy of body fat measurement: 0,1 %
- The range of water measurements in the body: 15-66,5 %
- Accuracy of water measurement in the body: 0,1 %
- Muscle tissue measurement range: 23-91 %
- Accuracy of muscle tissue measurement: 0,1 %
- Bone tissue measurement range: 1-6 %
- Accuracy of bone tissue measurements: 0,1 %
- Accuracy of the calculation of the recommended daily calorie intake (BMR): 1 kcal
- Power supply: 4x 1.5 V micro alkaline AAA batteries (LR03)
- Certification CE, RoHS
- Dimensions (w × h × h): 28 × 28 × 2.3 cm
- Weight (without batteries): 1465 g

PURPOSE OF USE

This product (personal scale) is used to measure the weight of the human body and with the help of a mobile device (phone) app, other biometric data and values such as metabolic age, body fat percentage, muscle and bone tissue in the body, etc.

PRINCIPLE OF MEASUREMENT

This scale uses bioimpedance analysis (BIA) to measure and analyse body parameters (percentage of bone, fat and muscle tissue, water in the body, etc.) based on the principle of a weak electrical impulse sent through the scale's foot sensors (electrodes) into the human body. The software embedded in the scale evaluates the measurements based on the difference in electrical impedance of the individual tissues and their relative contribution to the body.

Muscle contains about 90% water and adipose tissue 45% water. Therefore, the percentage sum can exceed 100% and for this reason the measurements cannot be added because both muscle and fat contain water.

SERVICE

Important notices:

- If the scale is placed on an uneven or soft surface (e.g. carpet), inaccuracies in weighing may occur (the scale may display higher or lower readings).
- We recommend weighing on the same pad/floor and at the same time of day, ideally first thing in the morning on an empty stomach - this will give comparable results.
- Wet, dirty feet or feet covered with cream can distort the results of body measurements or even make them impossible.
- Deviations in the measurement of body parameters may be greater in children under 10 years of age, persons over 70 years of age and professional athletes.
- Measurement of body parameters (except weight) is not recommended immediately after waking up, as the body is dehydrated after sleep. Similarly, we do not recommend measurements immediately after sports performance or during drastic weight loss diets.
- The measured values of the body parameters are only indicative and no medical conclusions can be drawn from them. If the results are significantly different from normal, we recommend consulting a nutritionist or doctor.
- The use of body parameter analysis mode (i.e. barefoot weighing on sensors) is not recommended:
 - pregnant women
 - people with a pacemaker or other medical implant
 - people with symptoms of oedema (swelling)
 - persons undergoing dialysis treatment

Inserting/replacing batteries

- 1) Open the A6 battery compartment cover at the bottom of the scale.
- 2) Insert the 4 AAA micro button cell batteries (LR03) included in the battery compartment. Make sure to observe the correct polarity, which is marked on the bottom of the battery compartment.
- 3) Close the battery compartment by snapping its cover back on.

Warning:

- Do not use rechargeable batteries to power the scale, their voltage may not be sufficient to power the scale.
- Do not combine old and new batteries due to their different residual capacities. A low battery error indication could be displayed prematurely and all batteries may need to be replaced.

Preparation for weighing

- 1) Place the scale on a flat, dry and hard surface.
- 2) The scale is equipped with the "Double Auto On" function, where the weighing can start immediately after you get on the scale. There is therefore no need to step on, switch on or otherwise activate the scale in advance.
- 3) The measured weight remains displayed on the display for about 25 seconds even after you get off the scale.

Switching weight units

This scale allows weighing in both metric and imperial systems. If necessary, switch the units using the switch at the bottom of the scale (**kg** → **lb** → **st lb**).

- The scale is factory preset to weigh in kilograms (the "kg" symbol appears next to the measured weight on the display).
- The weight unit you set remains stored in the scale's internal memory even after the scale is switched off and on again. The exception is when removing or replacing the battery, when the scale is switched back on, the unit for weighing the weight in grams is displayed (the g symbol appears on the display).


Weight weighing mode (without measuring body parameters)

In this mode, the scale weighs only the body weight.

- 1) Step onto the scale with your feet so that the entire surface of your feet touch the weighing surface. Do not stand on your toes or heels only, do not lean over and stand still until the resulting weight on the display is stable.
- 2) The resulting weight is indicated by a triple flash of its value on the display.
- 3) The resulting weight remains displayed for about 25 seconds, even after you have stepped off the scale. Then the display automatically switches off and after a further 20 seconds of inactivity the whole scale automatically switches off.

Weighing mode with body parameters measurement

In this mode, the scale weighs the weight and evaluates the body parameters at the same time, and it is necessary to stand on the scale with bare feet during weighing, otherwise the measurement cannot be made. The weighing surface is made of special glass, which allows the reading of the parameters of the human body through the feet even without the commonly used foot pads. Therefore, you must stand on the scale with the front of your foot at the top of the scale behind the dividing bar and the heel at the bottom of the scale below the dividing bar (see picture below).

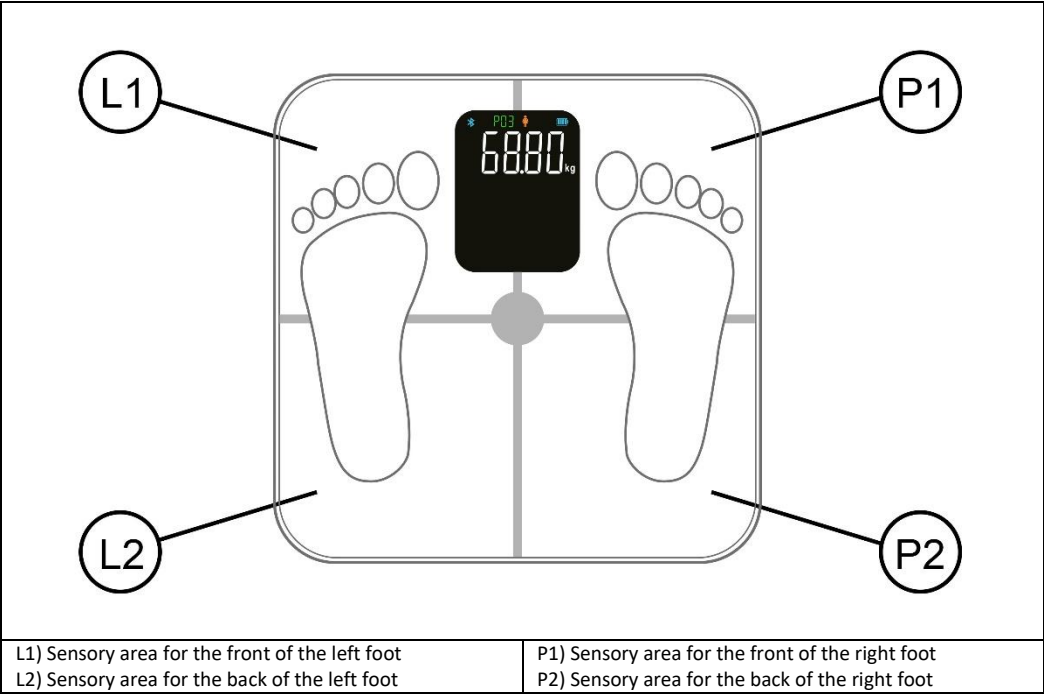- 1) Step onto the scale with bare feet so that the entire surface of the feet is touching the weighing surface. Do not stand on your toes or heels only, do not lean and stand still until the final weight is displayed.
- 2) The weighing progress is indicated by three "o" signs at the bottom of the display for each of the six parameters (see figure below).

- 3) The resulting weight and parameter values remain displayed for about 25 seconds. Then the display automatically switches off and after a further 20 seconds of inactivity the whole scale automatically switches off.
- 4) The values are further transferred to the application where they are analyzed and displayed.

MOVINGLIFE APP

This scale allows the transfer, analysis and display of all measured values to the MovingLife mobile app, which can be downloaded for free from the App Store (for Apple mobile phones) or Google Play (for Android mobile phones) by scanning the QR code below with your mobile phone. The values are transmitted wirelessly to the app via Bluetooth.

Warning:

The app requires Bluetooth version 4.2 support and, for mobile phones, Apple iOS 12.5 and above or Android 9.0 and above.

How to pair your mobile phone with a smart scale

- 1) On your mobile phone, turn on Bluetooth, connect it to a cellular and open the MovingLife app.
- 2) Press the weighing surface of the scale with your foot until the display reads "0.00 kg".
- 3) Within 25 seconds at the latest, select **I (Account) → My device → Add device (PLUS symbol)** on the mobile phone app. Select a **Salente BodyFit** personal scale and pair it.

Local account:

No registration is required, all data will only remain on your device. If you uninstall the app or start using another device, you will lose data.

Cloud account:

By switching the switch to the position, you can register to the Cloud and the data from the app will remain there even if you replace your device.

For the app to work properly and count health functions, you must fill in the basic user information that the app asks you during initial setup. The app can be linked to Apple Health or Google Fit, and synchronisation can be enabled in the app: **Me (Account) → Settings** and select Google Fit or Apple Health depending on your device.

Pairing failed

Please check the following items:

- 1) Check that the mobile phone operating system and Bluetooth version meet the compatibility conditions.
- 2) Stand on the scale and check that the display has lit up. If not, replace the batteries.
- 3) Make sure your mobile phone has Bluetooth capability.
- 4) If all the above items are OK, restart the app or restart your phone.

Options for switching off the scale

If the scale remains idle for more than 25 seconds, the display will automatically switch off. After 45 seconds of inactivity, the entire scale automatically shuts down.

CLEANING AND MAINTENANCE

- Never clean the scale under running water, rinse it or immerse it in water.
- Use a soft cloth dampened with lukewarm water to clean the outside of the scale. Do not use harsh cleaners, thinners or solvents or the scale surface may be damaged.

Disposal: If the information leaflet, packaging or product label contains a crossed-out bin symbol, this means that - in the European Union - all electrical and electronic products, batteries and accumulators marked in this way must be placed in separate waste bins at the end of their life. Do not dispose of these products in unsorted municipal waste.

DECLARATION OF CONFORMITY

Abacus Electric, s.r.o. hereby declares that this Salente BodyFit personal scale meets the requirements of the standards and regulations that are relevant to the type of equipment and at the same time declares that the product complies with Directive 2014/53/EU. The full text of the EU Declaration of Conformity is available at the following website:

<http://ftp.salente.cz/ce>

Importer / Manufacturer

Abacus Electric, s.r.o.

Planá 2, 370 01, Czech Republic

E-mail: helpdesk@salente.eu

Made in P.R.C.

Copyright © Abacus Electric, s.r.o.

www.salente.eu

All rights reserved.

Design and appearance are subject to change without notice.

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