
Nokia Thermo

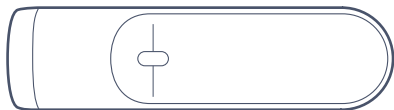QUICK INSTALLATION GUIDE

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Read this manual before use

EN | PL | CS | HU | SK

THANK YOU FOR CHOOSING NOKIA THERMO

DZIĘKUJEMY ZA WYBÓR NOKIA THERMO | DĚKUJEME, ŽE JSTE SI ZVOLILI VÝROBEK NOKIA THERMO
KÖSZÖNJÜK, HOGY A NOKIA THERMO KÉSZÜLÉKET VÁLASZTOTTA | ĎAKUJEME VÁM, ŽE STE SA ROZHODLI PRE PRODUKT NOKIA THERMO

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➔ nokia.com/health/support/thermo

This guide is available in additional languages at: www.nokia.com/health/guides

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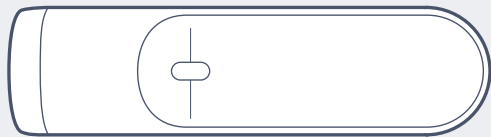
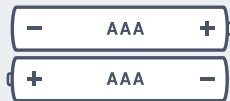
1

Visit thermoapp.nokia.com from your smartphone or tablet, and follow the instructions to install the Nokia Thermo app.

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
2

In the Thermo app, follow the instructions to install Thermo.

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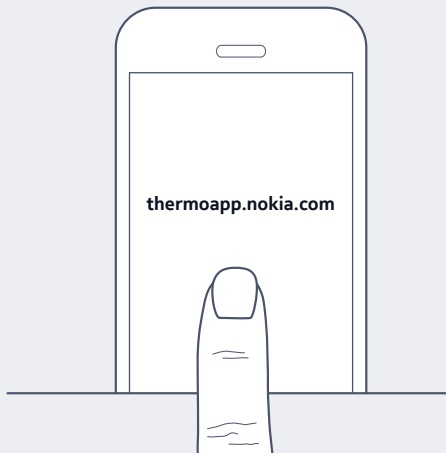1


Press the button to turn the thermometer on.

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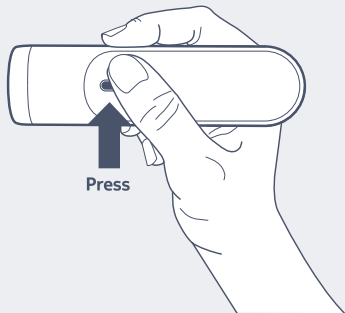Clear any sweat or hair from the forehead before taking a measurement.

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2

Starting from the middle of the forehead, slowly scan straight across to the hairline, up to 1cm above the skin.

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3

Thermo will vibrate twice when the measurement is complete.

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4


To assign the temperature to a user, scroll through names on the touch sensitive zone (1). Press the button to assign the reading (2).

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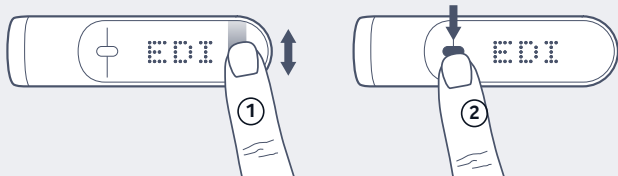5


A colored LED indicates the fever level according to the user's age.

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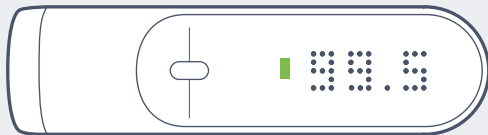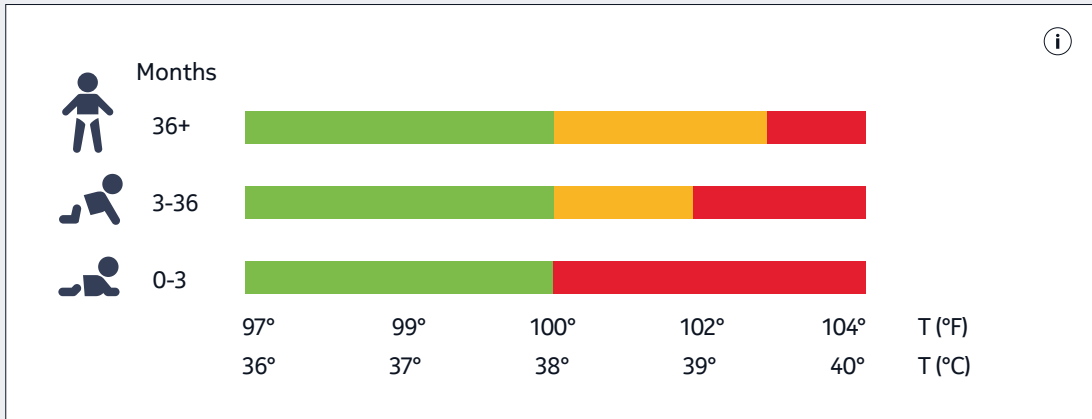
12

Person being measured must be indoors for 15 minutes prior to taking a reading. | Person being measured must have dry skin. If needed, dry with a towel prior to taking a reading. | Person being measured and the thermometer should be in the same ambient temperature for at least 15 minutes prior to taking a reading. | Person being measured should not have exercised within 15 minutes prior to taking a reading. | Don't take a measurement while or immediately after nursing a baby.

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
F Battery cover

PL Pokrywa baterii
CS Kryt baterií
HU Elem fedele
SK Kryt batérií

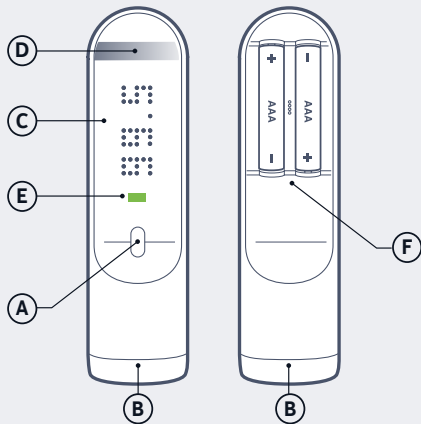US Nokia One (1) Year Limited Warranty - Nokia Thermo

Withings (Part of Nokia), 2 rue Maurice Hartmann, 92130 Issy-les-Moulineaux, France ("Nokia") warrants the Nokia-branded hardware product ("Nokia Product") against defects in materials and workmanship when used normally in accordance with Nokia's published guidelines for a period of ONE (1) YEAR from the date of original retail purchase by the end-user purchaser ("Warranty Period"). Nokia's published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Nokia does not warrant that the operation of the Nokia Product will be uninterrupted or error-free. Nokia is not responsible for damage arising from failure to follow instructions relating to the Nokia Product's use.

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UK Nokia One (1) Year commercial guarantee - Nokia Thermo

Withings (Part of Nokia), 2 rue Maurice Hartmann, 92130 Issy-les-Moulineaux, France ("Nokia") warrants the Nokia-branded hardware product ("Nokia Product") against defects in materials and workmanship when used normally in accordance with Nokia's published guidelines for a period of ONE (1) YEAR from the date of original retail purchase by the end-user purchaser ("Warranty Period"). Nokia's published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Nokia does not warrant that the operation of the Nokia Product will be uninterrupted or error-free. Nokia is not responsible for damage arising from failure to follow instructions relating to the Nokia Product's use.

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EN Withings (Part of Nokia) hereby declares that the device Nokia Thermo is in conformity with the essential requirements and other relevant requirements of Directive 2014/53/EU. The declaration of conformity can be found at: nokia.com/health/compliance

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NOKIA

Nokia Thermo | Smart temporal thermometer - Model: SCT01 - v1.0

© 2017 Nokia Technologies Ltd. All rights reserved.